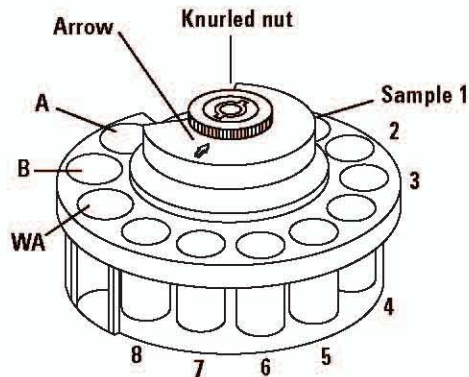


## Installing the 8-Vial Turret for the Automatic Liquid Sampler

### Agilent 7683B Automatic Liquid Sampler

You have just received the G2913-40010 8-vial turret used on the automatic liquid sampler (Figure 1). This 8-vial sample turret has two solvent positions (labeled **A** and **B**), and one waste position (labeled **WA**). You can control all necessary parameters using the 6890 series GC keyboard, the 6850 GC handheld controller, an Agilent data analysis system, or the G2911AA controller software. You cannot, however, use the G2614A or G2916A 100-vial sampler trays with this turret.



**Figure 1** 8-Vial turret

